
Reprocess EVIS EXERA™ III Duodenovideoscope TJF-Q190V with OER-Pro™ Endoscope Reprocessor

Quick Reference Guide

WARNING:
This cautionary note is only a summary of the steps necessary to properly reprocess TJF-Q190V. Be sure to follow the detailed steps outlined in the latest OPERATIONAL MANUAL that is included with your OER-Pro™ endoscope reprocessor.

Using CONNECTING TUBE MAJ-2358


1. Be sure to reprocess the TJF-Q190V with CONNECTING TUBE MAJ-2358 in OER-Pro endoscope reprocessors. (Figure 1)

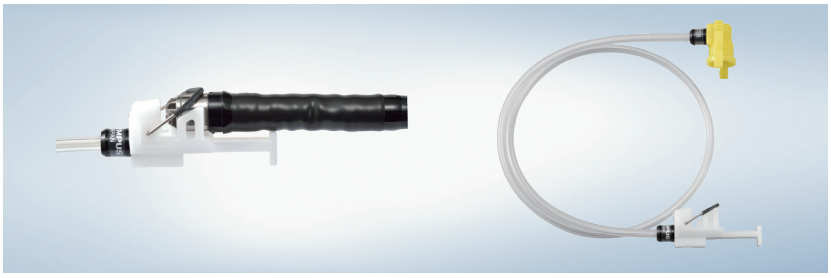


Reprocessing individually (NOT together with any other endoscope)


2. When reprocessing TJF-Q190V, be sure to reprocess individually. (Please do not reprocess TJF-Q190V in combination with any other endoscope.)

TJF-Q190V setup with MAJ-2358

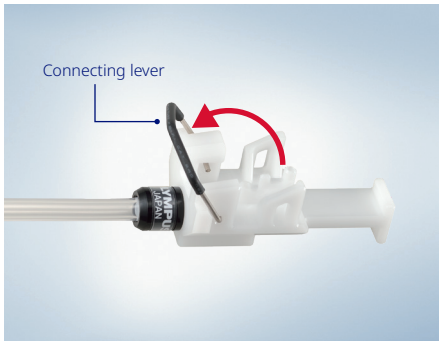

1. Raise the connecting lever

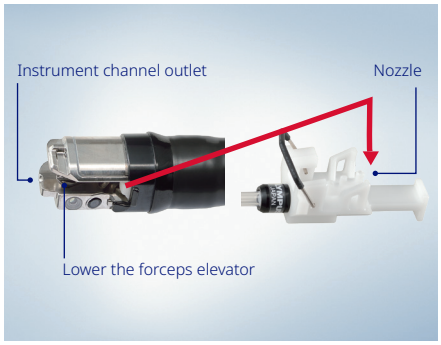

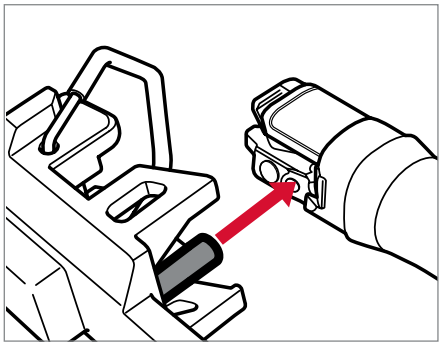


2. Insert the endoscope distal end by aligning the nozzle of the endoscope side connector into the instrument channel outlet of the endoscope


Reprocess EVIS EXERA™ III Duodenovideoscope TJF-Q190V with OER-Pro™ ENDOSCOPE REPROCESSOR

Quick Reference Guide


TJF-Q190V setup with MAJ-2358 (continued)

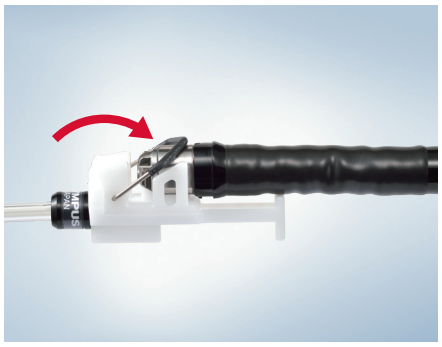


3. Lower the connecting lever to secure the distal end of the endoscope

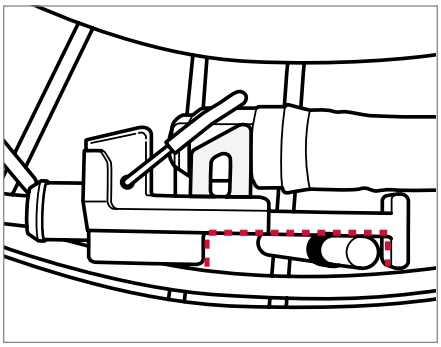



4. Set the distal end in the reprocessing basin while the groove of the endoscope side connector is attached to the pin on the proximal side of the retaining rack

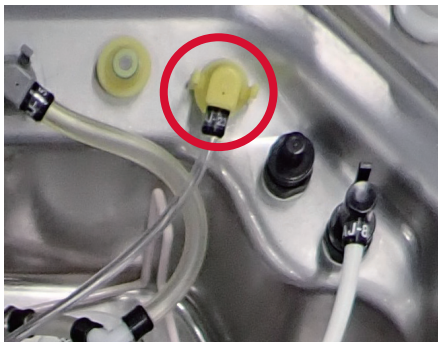

5. Connect the connecting tube with auxiliary water/elevator wire channel connector (yellow connector)

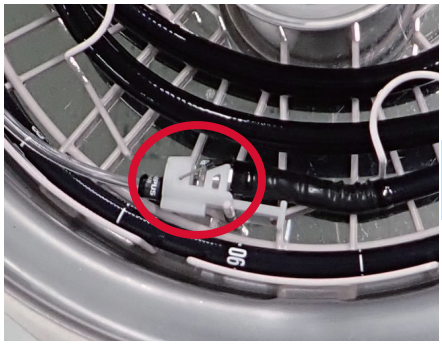

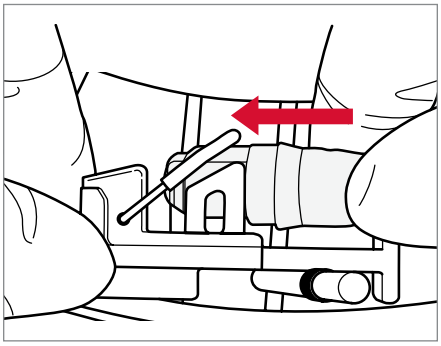

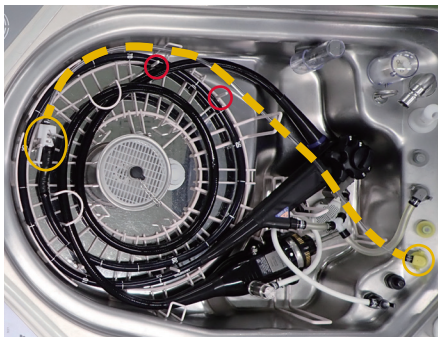

6. Gently push the endoscope distal end to ensure the nozzle of the endoscope side connector is in contact with the instrument channel outlet of the endoscope

* There are some versions of the OER-Pro endoscope reprocessor IFU that instruct to raise the forceps elevator to an angle of approximately 45° when reprocessing an endoscope which has a forceps elevator in an OER-Pro endoscope reprocessor. However, be sure to lower the forceps elevator when reprocessing the endoscope with the MAJ-2358 in OER-Pro endoscope reprocessors.

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.
Olympus is a registered trademark of Olympus Corporation, Olympus America Inc., and/or their affiliates. | Medical devices listed may not be available for sale in all countries.

OLYMPUS

OLYMPUS AMERICA INC.
3500 Corporate Parkway, PO BOX 610, Center Valley, PA 18034
www.medical.olympusamerica.com

For more information, contact your Olympus sales representative, or call 800-548-5515.

©2025 Olympus America Inc. All rights reserved.
Printed in USA OAICDS0925QRG57234